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INSTRUCTION MANUAL

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VZUSA

iCARE

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TABLE OF CONTENTS

- INTRODUCTION _____ p. 17
- SAFETY PRECAUTIONS _____ p. 17~18
- PRECAUTION BEFORE USE _____ p. 19
- PRECAUTION IN USE _____ p. 20
- PRECAUTION OF MAINTENANCE
AND INSPECTION _____ p. 21
- NAME OF PARTS _____ p. 22~23
- FUNCTIONAL DESCRIPTION _____ p. 24
- HOW TO CHARGE _____ p. 24
- HOW TO USE _____ p. 25
- TREATMENT PROCEDURE _____ p. 26~28
- PRODUCT SPECIFICATIONS _____ p. 29
- BATTERY RECYCLE _____ p. 29
- ASSURANCE PROVISION _____ p. 30
- WARRANTY CARD _____ Back cover

INTRODUCTION

- Thank you for purchasing the Vzusa iCARE.

In order to use the device properly and safely, it's necessary to read this instruction manual carefully and follow all safety precautions. Please include this instruction manual if you borrow or give to another person.

- This instruction manual contains the warranty card.

Confirm the information such as " The date of purchase" and " The purchase shop", and store it carefully,

SAFETY PRECAUTIONS

Be sure to read this instruction manual carefully before using, and always use properly.

To preserve the risk of injury, loss of life, and damage to any property, the following symbols are used. The marks and meanings of symbols are shown as below.

 DANGER	Failure to observe DANGER signs will result in hazard to personnel is high likely to result in serious injury or death.
 WARNING	Failure to observe WARNING signs could result in hazard to personnel possibly resulting in serious injury or death.
 CAUTION	Failure to observe CAUTION signs could result in injury to personnel and damage to the device and property. ※Property damage includes damage to house, household belongings, domestic animals and pets.

■ MEANINGS OF SYMBOLS

	This symbol indicates prohibition. ※It is indicated a special instruction inside or near a picture.
	This symbol indicates requirement that must be followed.
	This symbol indicates caution.

SAFETY PRECAUTIONS

WARNING

To preserve accidents or discomfort.

Prohibition

■ Do not use this product with any of the following medical electrical equipment.

Malfunction of medical electronic equipment causes accidents or discomfort.

- Biological implant such as pacemaker.
- Life support electronic equipment such as artificial heart lung.
- Wearing type medical electronic equipment such as electrocardiograph or drip infusion.

■ The following persons must not use the product.

- Persons with heart disease, abnormal blood pressure or malignant tumors.
- Persons with diseases with fever, infectious diseases or acute diseases.
- Persons who are receiving dental treatment.
- Persons with weakened immune systems or who use immune suppressants due to immunity disorders like AIDS etc.
- Persons with hemophilia, bleeding disorders, diabetes or facial neuralgia.
- Persons with paresthesia or dysalgesia.
- Persons with Keloid diathesis, facial melanosis, photosensitivity or pigmentation disorders.
- Persons with capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction.
- Persons who have allergies, atopic dermatitis, skin diseases, sensitive skin or weak skin.
- Persons who have epileptic stroke or discomfort due to flashing lights.
- Persons who have aftereffects of head or ophthalmologic aftereffects due to traffic accidents.
- Persons who have had eye operations (in the case of LASIK, consult with the doctor before use).
- Persons who are pregnant, on period or breast-feeding.
- Children or the person who cannot express their feelings in words.
- Persons under the influence of alcohol or with a hangover.
- Persons who take drugs that cause drowsiness.
- Persons who are extremely tired or do not feel well.
- Persons who have had skin problems such as spots or rashes due to artificial fibers.
- Persons who have or have had skin-inflammation due to cosmetics.
- Persons immediately after having chemical peeling at beauty salons etc.
- Persons who have had skin problems after using beauty equipment.

Prohibition

■ Do not use on the following parts.

- Except specific parts.
- Around heart, eyeballs, upper eyelids, lips, and mucous membrane such as inside the mouth.
- Private parts, anus, around Adam's apple, and ears.
- Parts with morbid spots, moles, and viral warts.
- Tinged and inflamed parts due to excessive sunburn.
- Parts with wounds, eczema, psoriasis, herpes, purulent acnes and/or pimples.
- Parts with pathological spots due to itching, flushing, physical stimulus.
- Parts where surgery was performed.
- Any parts immediately after hair removal or shaving has been performed.

Prohibition

※ If you are seeing a doctor, consult with the doctor before use.

PRECAUTION BEFORE USE

WARNING

- Do not use this product with other beauty equipment at the same time.
- Do not use this product for any purpose other than beauty and under conditions of malfunction or damage. This product is a treatment equipment specialized for use in beauty. Do not use for other purposes such as therapeutic purposes.
- Do not look at or directly view the LED light.
- Do not use continuously on the same parts.
Problems or symptoms such as skin damage, infection, inflammation, bleeding, pigmentation may be aggravated.
- Do not use serums that contain alcohol ingredients. It may cause skin problems.
- Do not use on pets.
- Do not use near flammable substances such as thinner, benzene, polish removers or sprays.
- Do not use with anything other than the dedicated adapter and USB cord. Do not connect the dedicated adapter and USB cord to other equipment. "The cause of conflagration, electrification, or breakdown."
- Do not damage the USB cord.
- ※ Do not damage, bend more than necessary, pull, twist, place heavy objects on top of, interpose. "Causes of fire, electric shock or damage."
- Do not adopt a star-burst connection. "Causes heat generation or fire."

Prohibition

Prohibition of use in a wet area

Prohibition of getting wet

Prohibition of use with wet hands

Unplug

Compulsion

- Do not use and keep in the bathroom or high humidity places.
- Do not leave in places where water is easily splashed "on a sink etc.".
- Do not immerse in water, pour or spill water (the causes of electric shock, fire or damage).
- Do not touch with wet hands. (the causes of electric shock or injury)
- Make sure to hold the end of the plug when putting in or pulling the adapter. Do not grab the cord to put in or pull the adapter.
- Do not use when the adapter is damaged or the insertion to any outlet is loose.
- Make sure to use only household outlets.
- Insert firmly up to the root of the adapter and the plug.
- Be sure to observe the contents of the instruction manual. Failure to do so may result in some problems.

PRECAUTION IN USE

WARNING

Unplug

- Stop using the product, turn off the power, and pull the adapter from the outlet in the following cases.
 - Earthquakes or lightning.
 - Smoke, abnormal heat, odor or distortion of the main body, the adapter or the code.

Prohibition

- During use, this product may slightly heat up but this is not a safety issue.
- Do not bring this product into contact with metals and any other electric appliances such as mobile phone, television.
- Charge completely before first time use this product.
- ※ To avoid turning off during use in case use no enough charge.
- Avoid to use in charge.

Compulsion

- Make sure to remove all metal objects, glasses and contact lenses during use.
- ※ Metal objects: Pierced earrings, necklaces, watches etc.
- Make sure to set to level 1 when using for the first time or until you get used to it.
- Do not press the head of product to the skin strongly and never place on the same part of the skin for longer than 3 seconds. This may cause skin troubles.
- During treatments around eyes, use this product with your eyes closed.
- Discontinue use immediately if abnormalities appear on the skin or in the body during use and consult with a doctor.

CAUTION

Precautions of the use of product.

Read the following sections carefully before use.

1

- Perform a special treatment to each part once a week.
- Perform a treatment twice a week in the first and second weeks. From the third week, perform a treatment about once a week on your skin conditions.
- Excessive treatment is not good for the skin.
- Give daily treatments to each part every day or every other day.
- Use this product according to your skin conditions.

2

- As to facial treatments, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin and apply the dedicated beauty essence or gel, then use this product.

3

- Skin conditions are different depending on your physical conditions.
- Use this product while observing your skin conditions, and discontinue use if abnormalities such as rashes, redness or itchiness appear on skin and consult with dermatologist.

PRECAUTION OF MAINTENANCE AND INSPECTION

WARNING

Prohibition of disassembly

- Do not disassemble, modify, or try to repair the device yourself. For repair, please inquire at dealer.

Prohibition

- Do not adhere dust and pin on AC adapter and plug. (Cause fire, electric shock or damage).
- Do not drop or hit this product.
- Do not wrap the code around the adapter when storing the device.
- Please keep out of reach of children and do not let them use this product.
- Do not store this product in the following places.
 - Places subject to direct sunlight, high humidity and dust.
 - Places subject to high temperature or temperature below freezing point or inside a car under the blazing sun.
 - Places subject to danger of ignition, rain or wind.
 - Places subject to chemicals, gas generation or excessive vibration.
 - Places being close to heating appliances, appliances generating heat and magnetic sources.

Compulsion

- How to clean
 - Wipe this product using a soft cloth dampened with water. If wiping with a hard cloth, the product may be damaged or this may cause deterioration in performance.
 - ※ Do not use solvents such as alkaline agents, alcohol, thinner or benzene, polish removers, polishing powders, hot water or chemicals.
 - ※ Do not leave this product with the beauty essence or gel on it after use.


INSPECTION

When this product is unused for longer than 6 months and is used again, make sure to inspect before use if there are any abnormalities (the lighting state of the LED lights and the presence of absence of cracks).

If used under faulty conditions, this may cause an accident.

NAME OF PARTS

Main Body

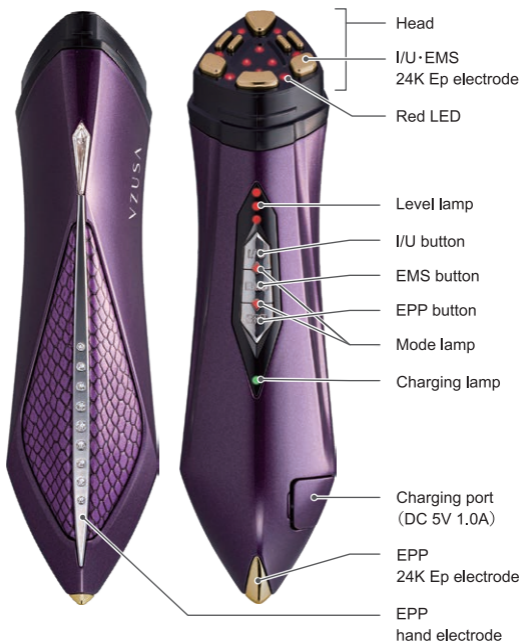Accessories

- AC adapter
- USB cable (Pin type)
- Cover
- Instruction manual with warranty card

NAME OF PARTS

24K Ep electrode

- The surface is surrounded with 24-karat gold which is more gentle on the skin.
- ※ This metal is really soft and easily scratched. Please handle it with sufficient care. If there are any small scratches, it will not any problem with quality of the device. It is safe to use.


Model equipped with NFC Tag

- To confirm Serial No., hold the NFC-compatible model on this product.

Compatible models

■ iOS : iPhone XS/XR or newer ■ Android : all

- ※ Models iPhone 5S/5C to iPhone X require a compatible app.
(Search for "NFC Reader" in the Apple store)
- ※ Not available with NFC-incompatible models.

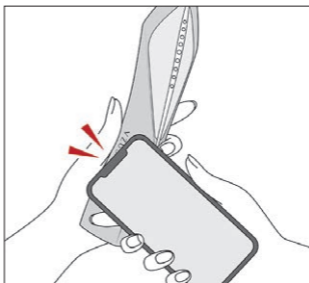FUNCTIONAL DESCRIPTION

Function

Button	Mode	Function
I/U	I/U mode (IPR / ULTRA PULSE)	IPR mode and UP mode output alternately. IPR (interval pulse rotation) = ULTRA PULSE outputs intermittently. ULTRA PULSE = Insert pulse / Mid-high frequency / EMS outputs simultaneously.
EMS	EMS mode	Output a pulse specialized in working muscles.
EPP	EPP mode	Give electrical stimulation to muscle and relax facial muscle.

LED

- Red LED will light up when starting the I/U mode or EMS mode. Red LED leads to elastic and lustrous skin.
- ※ Do not look directly at the LED light.
- ※ Keep eyes closed during treatment around eyes.

HOW TO CHARGE

- ※ Please fully charge the device when first use.
- ※ Do not use the device during charge, because it may accelerate battery degradation.
- ① Connect the AC adapter with USB code.
- ② Remove the cover of charging port, and insert USB cable into it.
- ③ Plug the AC adapter into a power outlet.
- ※ The lamp will light up when the battery starts charging.

■ Charge lamp (Green)

Low battery	Charging	Fully charged
Blinking	Lighting up	Lights out

HOW TO USE


- Make sure the battery is fully charged before use.

1 Mode

- Hold down the <I/U button> to start the I/U mode.
- Hold down the <EMS button> to start EMS mode.
- Hold down the <EPP button> to start the EPP mode.
- Level Lamp / Mode lamp will light up and the color of lamp differs depending on each mode.
 - RED : I/U mode ● Blue : EMS mode ● Green : EPP mode
- ※ Each mode starts at level 1.
- ※ I/U mode or EMS mode : Emit sound about once every 10 seconds since you start.
- ※ EPP mode : Emit sound when outputting.
- ※ To change the mode, long press the button of the mode on use to turn off the power and then reset again.

2 Level

- This product has 3 types of levels.
- Each level can be selected by pressing the button of the mode during use.
- ※ It is recommended to start using each mode with level 1 for the first time.

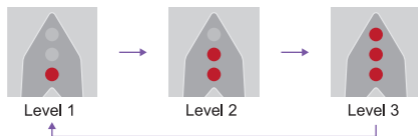3 Turn off the power


- Hold down the button of the mode during use to turn off the power after each treatment.
- ※ The device stops automatically in 10 minutes after starting up. (The sound can be heard once every second from before around 5 seconds of a stop.)

TREATMENT PROCEDURE

Before treatment

- 1 Remove accessories such as necklaces or rings and all metal objects that you are wearing.
- 2 Before performing a treatment, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin.

While performing a treatment

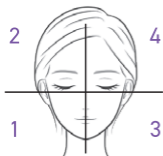■ When performing a treatment for full face, do not apply the beauty essence to the face all at once but divide it into four parts to apply as illustrated left and then perform a treatment.

EMS mode ▶ 10 minutes for daily treatment.

I/U mode ▶ 10 minutes for special treatment once a week

- 1 Apply the beauty essence on the treatment area and then perform a treatment.
- ※ During the treatment, apply the gel more if the treatment area becomes dry. When the treatment area becomes dry, it may cause irritation on your skin and the reduction of its effect.
 - ※ Learn more about treatment procedure, please see the page 27.

EPP mode ▶ 10 minutes for daily treatment.

- 1 Essence is not necessary, but if you care about the stimulation, you can spray some water onto the treatment areas and then perform a treatment.
- ※ Learn more about treatment procedure, please see the page 28.
 - ※ Do not use the device many times in a day cause it may hurt your skin.
 - ※ During the treatment, if the surface becomes rough, please check the device and clean up with soft cloth or tissues.
 - ※ Do not use the device if you are physically unwell.

After treatment


- 1 After treatment, gently spread the beauty essence and wipe away the gel left on the skin.
- 2 After treatment, wipe away the product with a soft cloth or tissues.

TREATMENT PROCEDURE

I/U mode

Special treatment : 10 minutes

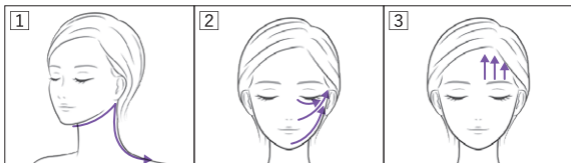
■ Apply beauty essence to the treatment area.

- 1 Move the device from the chin below to the rear side of the ear and then move down to the shoulder.
 - 2 Move the device from the end of the chin to the temple. Next, move the device from the nostril to the temple. Finally, move from inner corner of eye ,under eye, and to temple.
 - 3 Move the device from the eyebrows to the hairline.
 - Perform each procedure from 1 to 3 in the illustration for two or three times.
 - Move the device slowly and do not put the device on the same spot more than 3 seconds.
- ※ When the treatment area becomes dry, apply the beauty essence more.

EMS mode

Daily treatment : 10 minutes

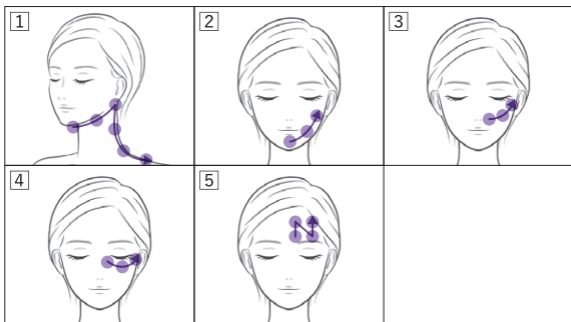■ Apply beauty essence to the treatment area.


- Use the device on each treatment point for approximately 3 seconds and then move to the next point.
 - Perform each procedure from 1 to 5 in the illustration for two or three times.
- ※ Keep the surface of all eight electrodes close to the treatment area.
- ※ When the treatment area becomes dry, apply the beauty essence more.
- ※ The procedure from illustration is a guideline only.

TREATMENT PROCEDURE

EPP mode

Daily treatment : 10 minutes

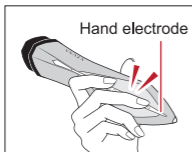- Hold the hand electrode of the main body with bare hands.
- ※ EPP mode will not be able to output properly if the hand electrode is not held by bare hands.

- Keep the device on the skin and move to the next spot in every 5 seconds following the directions of the illustration.
- Perform each procedure from 1 to 3 in the illustration two or three times.
- ※ Essence is not necessary, but if you care about the stimulation, you can spray some water onto the treatment areas and then perform a treatment.
- ※ The device emits "beep" sound while EPP mode is outputting.
- ※ It is not matter to feel electricity to a hand touching earth plate.
- ※ When the treatment area becomes dry, apply the beauty essence more.

WARNING

Do not use around mucosa, upper eyelids, eyes, ears, and Adam's apple.

PRODUCT SPECIFICATIONS

Product name	the Vzusa iCARE
Power rating of AC adapter	Input 100-240V ~ 50/60Hz Output 5V = 1.0A
Power rating of Main body	Input 5V = 1.0A
Power consumption	Maximum 6W
Size	Main body+Cover : Approx. H190×W50×D46mm
Weight	Main body+Cover : Approx. 115g

BATTERY RECYCLE

Make sure to take out the lithium ion battery when disposing the device.

Make sure to consume the power in the lithium ion battery before taking out the the lithium ion battery.

Do not discard the the lithium ion battery. Please follow your local government rules to dispose in the specified battery recycling place.

(Battery removal procedure:Unclench the shell of the device with a tool, then take out the battery. Be careful to avoid the the battery's positive and negative short-circuiting.)

Li-ion

ASSURANCE PROVISION

■ This warranty promises free repair if the conditions written below are fulfilled. If the product is damaged within the warranty period from the date of purchase, please contact the purchase shop or our customer center by submitting a copy of this warranty.

※ In the case that the serial number, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card, the warranty is not valid so please make sure to check that you have not missed out any required information in advance.

※ Accessories (AC adapter, USB cable, Cover and instruction manual with warranty) are not covered by the warranty.

※ In the case that the serial No. stickers on the main body and the warranty card are removed, the product is not covered by the warranty.

■ **The repairing cost will be charged in the following cases even within the warranty period.**

A. In the case of the removal of the serial No. stickers on the warranty card.

B. In the case the serial No, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card.

C. Malfunction or damage due to improper handling (e.g. dropping and submerging in water) caused during transportation or moving by the user.

D. Malfunction or damage due to fire, other natural disasters or abnormal voltage.

E. Damage due to other devices connected to this product.

F. Damage due to the use of an AC adapter other than the be dedicated AC adapter included with this product.

G. In the case of repair, adjustment or modification by anyone other than us.

H. In the case of not submitting of the warranty card.

■ **This warranty card is valid only in Japan.**

※ In the case of purchasing this product overseas, the warranty covers the product only when the shop where user purchased the product has an authorized distributor certificate issued by manufacturer.

■ **Shipping fee**

a. Our company pays the round delivery costs during warranty period. "Be effective only in Japan"

b. The round delivery cost will be charged to customer after the warranty period.

■ **This warranty card is not reissued. Please keep it carefully.**

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