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INSTRUCTION MANUAL

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# VZUSA

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# TABLE OF CONTENTS

- INTRODUCTION \_\_\_\_\_ p. 18
- SAFETY PRECAUTIONS \_\_\_\_\_ p. 18~19
- PRECAUTION BEFORE USE \_\_\_\_\_ p. 20
- PRECAUTION IN USE \_\_\_\_\_ p. 21
- PRECAUTION OF MAINTENANCE  
AND INSPECTION \_\_\_\_\_ p. 22
- NAME OF PARTS \_\_\_\_\_ p. 23~24
- EXPLANATION OF FUNCTIONS \_\_\_\_\_ p. 25
- HOW TO USE \_\_\_\_\_ p. 26
- TREATMENT PROCEDURE \_\_\_\_\_ p. 27
- TREATMENT PROCEDURE (Face) \_\_\_\_\_ p. 28
- TREATMENT PROCEDURE (Body) \_\_\_\_\_ p. 29~30
- PRODUCT SPECIFICATIONS \_\_\_\_\_ p. 31
- ASSURANCE PROVISION \_\_\_\_\_ p. 32
- MEMO \_\_\_\_\_ p. 50
- WARRANTY CARD \_\_\_\_\_ Back cover

# INTRODUCTION

- Thank you for purchasing the Vzusa.

In order to use the device properly and safely, it's necessary to read this instruction manual carefully and follow all safety precautions. Please include this instruction manual if you borrow or give to another person.

- This instruction manual contains the warranty card.

Confirm the information such as " The date of purchase" and " The purchase shop", and store it carefully.

# SAFETY PRECAUTIONS

**Be sure to read this instruction manual carefully before using, and always use properly.**

To preserve the risk of injury, loss of life, and damage to any property, the following symbols are used. The marks and meanings of symbols are shown as below.

 <b>DANGER</b>	Failure to observe DANGER signs will result in hazard to personnel is high likely to result in serious injury or death.
 <b>WARNING</b>	Failure to observe WARNING signs could result in hazard to personnel possibly resulting in serious injury or death.
 <b>CAUTION</b>	Failure to observe CAUTION signs could result in injury to personnel and damage to the device and property. ※Property damage includes damage to house, household belongings, domestic animals and pets.

## ■ MEANINGS OF SYMBOLS

	This symbol indicates prohibition. ※It is indicated a special instruction inside or near a picture.
	This symbol indicates requirement that must be followed.
	This symbol indicates caution.

# SAFETY PRECAUTIONS

## WARNING

### To preserve accidents or discomfort.

Prohibition

■ **Do not use this product with any of the following medical electrical equipment.**

- Malfunction of medical electronic equipment causes accidents or discomfort.
- Biological implant such as pacemaker.
- Life support electronic equipment such as artificial heart lung.
- Wearing type medical electronic equipment such as electrocardiograph or drip infusion.

Prohibition

■ **The following persons must not use the product.**

- Persons with heart disease, abnormal blood pressure or malignant tumors.
- Persons with diseases with fever, infectious diseases or acute diseases.
- Persons who are receiving dental treatment.
- Persons with weakened immune systems or who use immune suppressants due to immunity disorders like AIDS etc.
- Persons with hemophilia, bleeding disorders, diabetes or facial neuralgia.
- Persons with paresthesia or dysalgesia.
- Persons with Keloid diathesis, facial melanosis, photosensitivity or pigmentation disorders.
- Persons with capillary dilatation due to long-term steroidal hormone treatments or liver dysfunction.
- Persons who have allergies, atopic dermatitis, skin diseases, sensitive skin or weak skin.
- Persons who have epileptic stroke or discomfort due to flashing lights.
- Persons who have aftereffects of head or ophthalmologic aftereffects due to traffic accidents.
- Persons who have had eye operations (in the case of LASIK, consult with the doctor before use).
- Persons who are pregnant, on period or breast-feeding.
- Children or the person who cannot express their feelings in words.
- Persons under the influence of alcohol or with a hangover.
- Persons who take drugs that cause drowsiness.
- Persons who are extremely tired or do not feel well.
- Persons who have had skin problems such as spots or rashes due to artificial fibers.
- Persons who have or have had skin-inflammation due to cosmetics.
- Persons immediately after having chemical peeling at beauty salons etc.
- Persons who have had skin problems after using beauty equipment.

Prohibition

■ **Do not use on the following parts.**

- Except specific parts.
- Around heart, eyeballs, upper eyelids, lips, and mucous membrane such as inside the mouth.
- Private parts, anus, around the Adam's apple, and ears.
- Parts with morbid spots, moles, and viral warts.
- Tinged and inflamed parts due to excessive sunburn.
- Parts with wounds, eczema, psoriasis, herpes, purulent acnes and/or pimples.
- Parts with pathological spots due to itching, flushing, physical stimulus.
- Parts where surgery was performed.
- Any parts immediately after hair removal or shaving has been performed.

※ If you are seeing a doctor, consult with the doctor before use.

# PRECAUTION BEFORE USE

## WARNING

- Do not use this product with other beauty equipment at the same time.
- Do not use this product for any purpose other than beauty and under conditions of malfunction or damage. This product is a treatment equipment specialized for use in beauty. Do not use for other purposes such as therapeutic purposes.
- Do not look at or directly view the LED light.
- Do not use continuously on the same parts.  
Problems or symptoms such as skin damage, infection, inflammation, bleeding, pigmentation may be aggravated.
- Do not use serums that contain alcohol ingredients. It may cause skin problems.
- Do not use on pets.
- Do not use near flammable substances such as thinner, benzene, polish removers or sprays.
- Do not use with anything other than the dedicated adapter and USB cord. Do not connect the dedicated adapter and USB cord to other equipment. "The cause of conflagration, electrification, or breakdown."
- Do not damage the USB cord.
- ※ Do not damage, bend more than necessary, pull, twist, place heavy objects on top of, interpose. "Causes of fire, electric shock or damage."
- Do not adopt a star-burst connection. "Causes heat generation or fire."

Prohibition

Prohibition of use in a wet area

Prohibition of getting wet

Prohibition of use with wet hands

Unplug

Compulsion

- Do not use and keep in the bathroom or high humidity places.
- Do not leave in places where water is easily splashed "on a sink etc."
- Do not immerse in water, pour or spill water (the causes of electric shock, fire or damage).
- Do not touch with wet hands. (the causes of electric shock or injury)
- Make sure to hold the end of the plug when putting in or pulling the adapter. Do not grab the cord to put in or pull the adapter.
- Do not use when the adapter is damaged or the insertion to any outlet is loose.
- Make sure to use only household outlets.
- Insert firmly up to the root of the adapter and the plug.
- Be sure to observe the contents of the instruction manual. Failure to do so may result in some problems.

# PRECAUTION IN USE

## WARNING

Unplug

- Stop using the product, turn off the power, and pull the adapter from the outlet in the following cases.
  - Earthquakes or lightning.
  - Smoke, abnormal heat, odor or distortion of the main body, the adapter or the cord.

Prohibition

- During use, this product may slightly heat up but this is not a safety issue.
- Do not bring this product into contact with metals and any other electric appliances such as mobile phone, television.

Compulsion

- Make sure to remove all metal objects, glasses and contact lenses during use.
- ※ Metal objects: Pierced earrings, necklaces, watches etc.
- Make sure to set to level 1 when using for the first time or until you get used to it.
- Do not press the head of product to the skin strongly and never place on the same part of the skin for longer than 3 seconds. This may cause skin troubles.
- During treatments around eyes, use this product with your eyes closed.
- Discontinue use immediately if abnormalities appear on the skin or in the body during use and consult with a doctor.

## CAUTION

### Precautions of the use of product.

Read the following sections carefully before use.

1

- Perform a special treatment to each part once a week.
- Perform a treatment twice a week in the first and second weeks. From the third week, perform a treatment about once a week on your skin conditions.
- Excessive treatment is not good for the skin.
- Give daily treatments to each part every day or every other day.
- Use this product according to your skin conditions.

2

- As to facial treatments, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin and apply the dedicated beauty essence or gel, then use this product.
- Body treatment should be applied after bathing or in other clean skin conditions.

3

- Skin conditions are different depending on your physical conditions.
- Use this product while observing your skin conditions, and discontinue use if abnormalities such as rashes, redness or itchiness appear on skin and consult with dermatologist.

# PRECAUTION OF MAINTENANCE AND INSPECTION

## WARNING

Prohibition of disassembly

- Do not disassemble, modify, or try to repair the device yourself. For repair, please inquire at dealer.

Unplug

- Unplug the adapter from the outlet when cleaning the main body. (Failure to comply may cause fire or accident.)
- Make sure to turn off the power, and disconnect the AC adapter with the outlet when not using. (Failure to comply may cause fire or damage.)

Prohibition

- Do not adhere dust and pin on AC adapter and plug. (Cause fire, electric shock or damage.)
- Do not drop or hit this product.
- Do not wrap the code around the adapter when storing the device.
- Please keep out of reach of children and do not let them use this product.
- Do not store this product in the following places.
  - Places subject to direct sunlight, high humidity and dust.
  - Places subject to high temperature or temperature below freezing point or inside a car under the blazing sun.
  - Places subject to danger of ignition, rain or wind.
  - Places subject to chemicals, gas generation or excessive vibration.
  - Places being close to heating appliances, appliances generating heat and magnetic sources.

Compulsion

- How to clean  
Wipe this product using a soft cloth dampened with water. If wiping with a hard cloth, the product may be damaged or this may cause deterioration in performance.
  - ※ Do not use solvents such as alkaline agents, alcohol, thinner or benzine, polish removers, polishing powders, hot water or chemicals.
  - ※ Do not leave the device with the essence or gel on it after use.


## INSPECTION

When this product is unused for longer than 6 months and is used again, make sure to inspect before use if there are any abnormalities (the lighting state of the LED lights and the presence of absence of cracks).

If used under faulty conditions, this may cause an accident.

# NAME OF PARTS

## Main Body
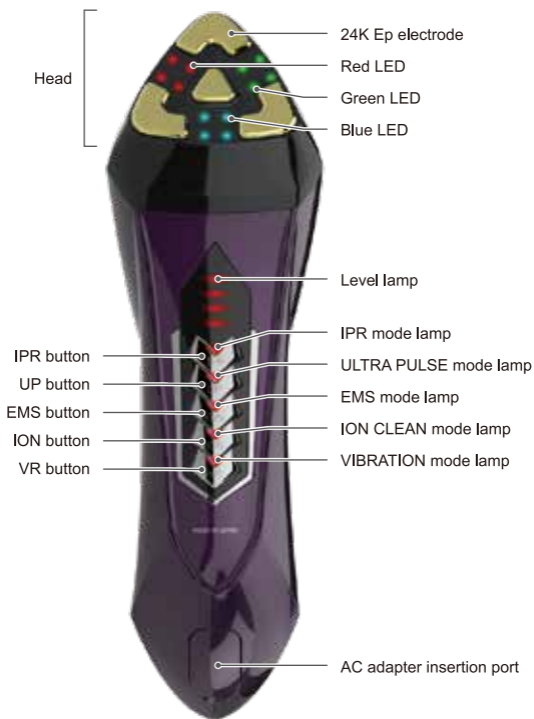
## Accessories

- AC adapter
- Stand
- Instruction manual with warranty card

# NAME OF PARTS

## 24KEp electrode

- The surface is surrounded with 24-karat gold which is more gentle on the skin.
- ※ This metal is really soft and easily scratched. Please handle it with sufficient care. If there are any small scratches, it is not any problem with quality of the device. It is safe to use.

## Model equipped with NFC Tag

- To confirm Serial No., hold the NFC-compatible model on this product.

### Compatible models

- iOS: iPhone XS/XR or newer
- Android: all

- ※ Models iPhone 5S/5C to iPhone X require a compatible app. (Search for "NFC Reader" in the Apple store)
- ※ Not available with NFC-incompatible models.

# EXPLANATION OF FUNCTIONS

## Function

Button	Mode	Function
IPR	IPR mode	Output Ultra Pulse Intermittently on skin and rotate quickly. IPR=Interval Pulse Rotation.
UP	ULTRA PULSE mode	Output insert pulse / EMS / mid-high frequency simultaneously.
EMS	EMS mode	Output a pulse specialized in working muscles.
ION	ION CLEAN mode	Absorb and eliminate impurities from skin with ionic power.
VR	VIBRATION mode	Mitigate electric stimulation with the vibration.

## LED

■ Three colors of LED lights with different functions can provide overall treatment.

The blue light leads to the healthy skin. The green light leads to the translucent and impressive skin. The red light leads to the elastic and lustrous skin.

- ※ Do not look directly at the LED light.
- ※ Keep your eyes closed during the treatment around eyes.

# HOW TO USE

※ When using this device for the first time, remove the protective sticker (the transparent sticker) from the electrode.

① Open the cover of the AC adapter socket of the main body and plug the AC adapter into the main body and the outlet.

※ Perform a treatment in the any mode of IPR mode, ULTRA PULSE mode, EMS mode, and ION CLEAN mode.

② Mode selection

- Press the 〈IPR button〉 to select the IPR mode.
- Press the 〈UP button〉 to select the ULTRA PULSE mode.
- Press the 〈EMS button〉 to select the EMS mode.
- Press the 〈ION button〉 to select the ION CLEAN mode.

※ The mode lamp of the selected mode is lit up in red.

※ Press the button and start a treatment with level 1.

※ A sound is made once every 10 seconds.


③ Level selection

- This product is 4 types of levels.

Choose level by pressing the button of the mode that selected.

※ When using this device for first time, start with level 1, choose the level that you like and then perform a treatment with that level.

※ When reducing the level or changing the mode, press the button of the mode on use to turn off the power. Then, reset the level or select the mode.

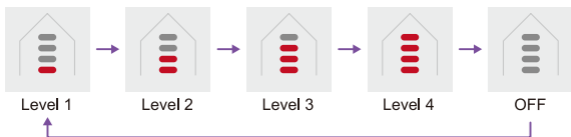④ VIBRATION selection

- Press the 〈VR button〉 to use VIBRATION mode as you like.
- VIBRATION mode has two patterns. The order of switching is Short ▶ Long ▶ OFF.

※ Press the 〈VR button〉 and only VIBRATION mode can be used.

⑤ Power OFF

- With the activation mode at level 4, press the button again to turn off the power.

※ The unit will automatically stop 10 minutes after startup. (A sound will be emitted 5 seconds before the shutdown.)

⑥ Disconnect the AC adapter from the main unit and outlet.


※ Close the cover of the AC adapter socket of the main unit after use.

# TREATMENT PROCEDURE

## Before treatment

- 1 Remove accessories such as necklaces or rings and all metal objects that you are wearing.
- 2 Before performing a treatment, remove makeup and wash your face to remove impurities such as makeup cosmetics and sebum on your skin.  
Body treatment should be applied after bathing or in other clean skin conditions.

## During treatment

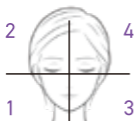- When performing a treatment for full face, do not apply the beauty essence to the face all at once but divide it into four parts to apply as illustrated left and then perform a treatment.

## Special treatment once a week

- **ION CLEAN mode, IPR mode and ULTRA PULSE mode = 15 minutes in total**
- 1 Spray the dedicated water on the part you want to treat, and perform a treatment with the ION CLEAN mode.
- 2 Wash your face or wipe off your skin entirely with a cotton pad.
- 3 Apply the beauty essence on the treatment parts, and perform a treatment with either the IPR mode or the ULTRA PULSE mode.
- ※ It is effective to perform 5-minute treatment with IPR mode and perform 5-minute treatment with UP mode.
- ※ Never place on the same part of the skin for longer than 3 seconds. This may cause skin problems.

## Daily treatment

- **EMS mode = 10 minutes**
- 1 Apply the beauty essence on the treatment parts.
- ※ Excessive treatment burdens the skin. Never use this product several times a day.
- ※ During treatment, apply more dedicated water and beauty essence if you feel dryness. When the application quantity is too little you may feel an irritation on the skin or the effect of the treatment may not be exhibited.
- ※ Place all four elements closely to the part of the skin to be treated.
- ※ During use, in the case that the surface of the product's head becomes rough by the adhesion of dirt etc., wipe the surface of the head using a soft cloth, tissue, etc.
- ※ Do not use this product when you do not feel well.

## After treatment

- 1 After treatment, gently spread the beauty essence left on the skin.
- 2 After treatment, wipe away the product with a soft cloth or tissues.


# TREATMENT PROCEDURE (Face)

## IPR mode / ULTRA PULSE mode / ION CLEAN mode

- 1 Move along the face line from below the jaw to the rear side of the ear. Then, move from the parotid gland to the shoulder.
  - 2 Move along the face line from the end of the jaw to the temple. Next, move along the cheek from the nostril to the temple. Finally, move from inner corner of eye, under eye, and to temple.
  - 3 Move from the bottom to the top while lifting the eyebrows upward to the hairline.
- Perform each procedure 1 to 3 two or three times.
  - Move constantly and slowly and do not stay on the same part for any length of time.
- ※ When the dedicated water or beauty essence becomes dry, apply more.

## EMS mode

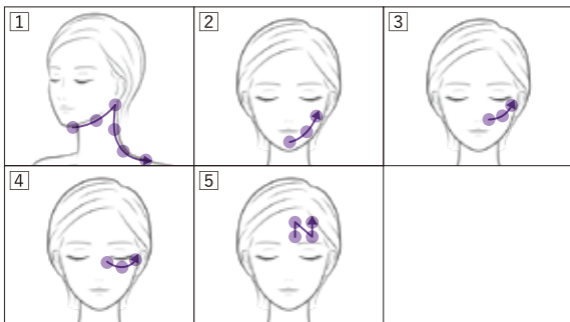- Perform each procedure 1 to 5 two or three times.
  - Place this product to the skin in order for approx. 3 seconds as shown in the illustration.
- ※ The treatment procedure in the illustration is reference.
- ※ When the beauty essence beauty essence becomes dry, add more.


### ⚠ WARNING

Do not use around mucosa, upper eyelids, eyeballs, ears, and Adam's apple.

# TREATMENT PROCEDURE (Body)

## Upper arm


- Special treatment: each arm 15 minutes
- Daily treatment: each arm 15 minutes

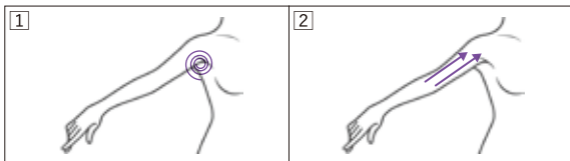- 1 Move like drawing small circles around the root of arm two or three rounds.
- 2 Move straight along the upper arm from above the elbow to the root of arm.

## Waist

- Special treatment: 15 minutes
- Daily treatment: 15 minutes

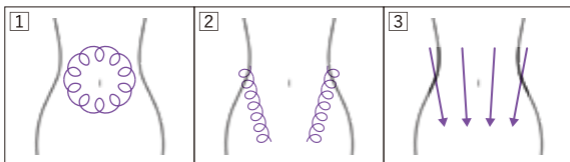- 1 Move like drawing small circles around the navel once round.
- 2 Move like drawing small circles from the side of the waist to the inguinal node as if making the curve of the waist.
- 3 Move straight as if flowing from the top to the bottom.

## Hip

- Special treatment: 15 minutes
- Daily treatment: 15 minutes


- 1 Move like drawing small circles from the boundary between the hipline and the thighs to below the waist.
- 2 Move in zigzags around the boundary between the hipline and the thighs.
- 3 Move straight from the boundary between the hipline and the thighs to below the waist.

# TREATMENT PROCEDURE (Body)

## Leg (front and back)


- Special treatment: each calf 15 minutes
- Daily treatment: each calf 15 minutes

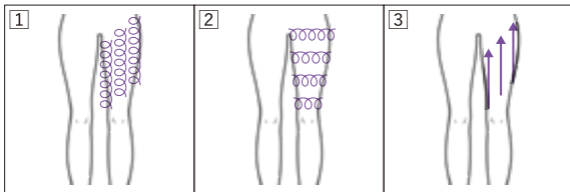- 1 Move sideways like drawing small circles from above the knee (the back of the knee) to the root of the thighs (the hipline).
  - 2 Move linearly like drawing small circles from above the knee (the back of the knee) to the root of the thighs (the hipline).
  - 3 Move straight from above the knee (the back of the knee) to the root of the thighs (the hipline).
- ※ Perform in the same way for the front and the back of the legs.

## Calf

- Special treatment: each calf 15minutes
- Daily treatment: each calf 15 minutes

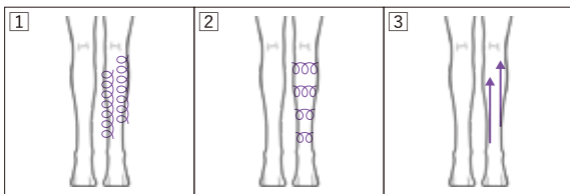- 1 Move like drawing small circles from the ankle to the back of the knee.
- 2 Move sideways like drawing small circles from the ankle to the back of the knee.
- 3 Move linearly from the ankle to the back of the knee.

# PRODUCT SPECIFICATIONS

Product name	the Vzusa
Power rating of AC adapter	Input 100-240V ~ 50/60Hz Output 12V = 1.5A
Power rating of Main body	Input 12V = 1.5A
Power consumption	Maximum 12W
Size	Main body: Approx. H208 × W63.4 × D60mm
Weight	Main body: Approx. 200g

# ASSURANCE PROVISION

■ This warranty promises free repair if the conditions written below are fulfilled. If the product is damaged within the warranty period from the date of purchase, please contact the purchase shop or our customer center by submitting a copy of this warranty.

- ※ In the case that the serial number, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card, the warranty is not valid so please make sure to check that you have not missed out any required information in advance.
- ※ Attachments (AC adapter, stand, instruction manual with warranty) are not covered by the warranty.
- ※ In the case that the serial No. stickers on the main body and the warranty card are removed, the product is not covered by the warranty.

■ **The repairing cost will be charged in the following cases even within the warranty period.**

- A. In the case of the removal of the serial No. stickers of the main body and the warranty card.
- B. In the case the serial No, date of purchase, name of dealer, customer name, address and phone number are not filled out on the warranty card.
- C. Malfunction or damage due to improper handling (e.g. dropping and submerging in water) caused during transportation or moving by the user.
- D. Malfunction or damage due to fire, other natural disasters or abnormal voltage.
- E. Damage due to other devices connected to this product.
- F. Damage due to the use of an AC adapter other than the be dedicated AC adapter included with this product.
- G. In the case of repair, adjustment or modification by anyone other than us.
- H. In the case of not submitting of the warranty card.
- I. Damage due to a beauty essence by using a beauty essence other than the dedicated beauty essence recommended by us.

■ **This warranty card is valid only in Japan.**

- ※ In the case of purchasing this product overseas, the warranty covers the product only when the shop where user purchased the product has an authorized distributor certificate issued by manufacturer.

■ **Shipping fee**

- a. Our company pays the round delivery costs during warranty period. "Be effective only in Japan"
- b. The round delivery cost will be charged to customer after the warranty period.

■ **This warranty card is not reissued. Please keep it carefully.**

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